



TE II

The **TE II** (Total Experience II Workstation) integrates the worksurface, standing area, and traffic lanes, using every inch of available floor space to maximum benefit. The unique design of the TE II eliminates the need for students to work back to back, creating a spacious, safe working environment.

Features


- Fire-retardant understructure of molded fiberglass bolted to a 10-gauge steel support pedestal
- Vandal-resistant Unicast CW/G fixtures are standard
- Raised turrets allow extra height
- Top has marine edge to control spills
- Several models are available
- Storage cubicle options are available
- Electrical options are available

sales@sheldonlabs.com



800.531.7604

TE II

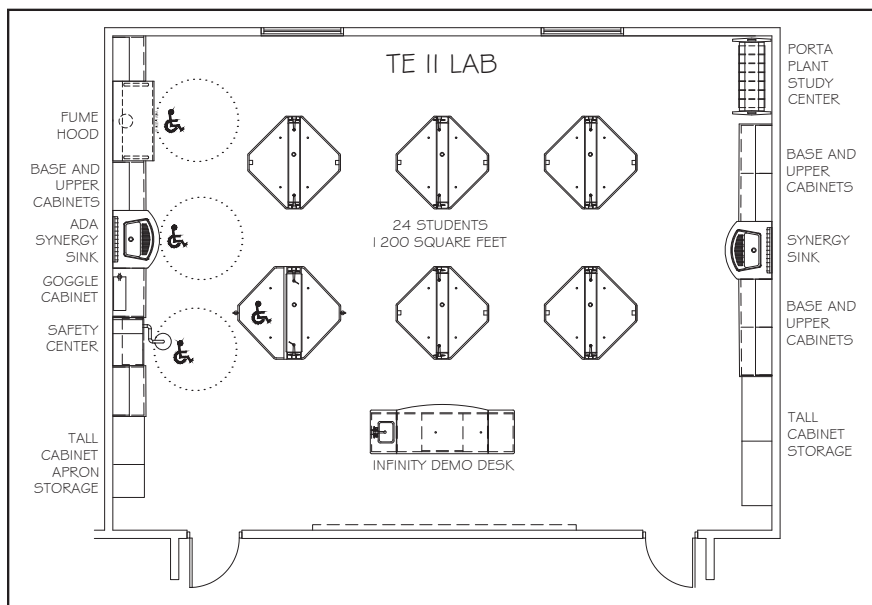

TE II Layout

The TE II table lab layout integrates the worksurface and traffic path, eliminating the need for students to work back-to-back, creating a spacious working environment. The TE II is an excellent choice for chemistry, physics, biology, and other science disciplines.

Worksurface

Epoxy resin with marine perimeter edge to contain spillage.

Sink Basin

Available in trough or round style (epoxy resin). Also available without sink.

Services

Unicast fixtures (cold water & double gas) are standard. Unimix fixtures (cold water, hot water and gas) are optional. Internal piping, wiring, conduit and drain lines are not included.

Storage

Available as mobile or attached units.

Upright Rods

Upright rod sockets included in top. Upright rods are optional.

Support Understructure

Molded one-piece fiberglass reinforced polyester understructure.

Steel Pedestal Base

Height adjustable at time of installation. Includes removable steel panels for access to services. Standard color is black.

Accessibility

Worksurface available with split top with one surface at ADA height (max. 34" high). Includes unicast fixture with an ADA wristblade handle.

Connectivity

Two (2) GFI duplex receptacles included on support understructure.


Table Options

Over 10 table models available; refer to TE II Matrix for more details.

Color Options-Epoxy

Available in Black Onyx, Lunar White, Gray, Tan, and Dark Khaki.

Color Options-Fiberglass*

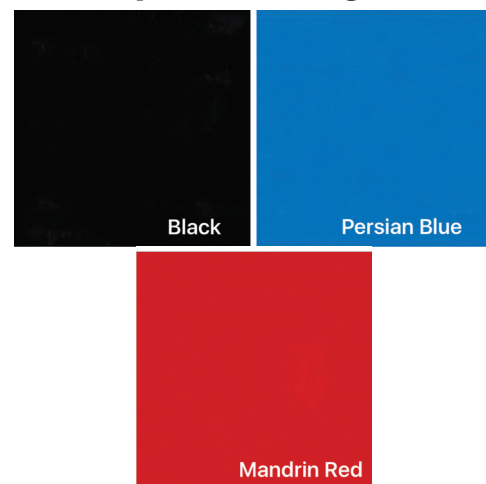

***If color match is critical, samples are available upon request.**

sales@sheldonlabs.com



800.531.7604